

EXPERIMENT - 1

Object :- To prepare and standardize the solution of NaOH (N/10) against standard solution of Oxalic acid (N/5).

APPARATUS/CHEMICALS : Burette, Pipette, conical flask, NaOH solution, phenolphthalein indicator, oxalic acid (N/5).

THEORY :

Standard solution is one in which amount of a substance is present in definite volume of the solution, whose concentration is known called as standard solution.

Types of solution

- 1) Primary standard solution
- 2) Secondary standard solution

PRIMARY STANDARD SOLUTION :-

The substance whose standard solution is prepared by dissolving directly its definite amount and then dissolving in definite volume of solvent or solution is



Name of Experiment..... Date.....

Experiment No..... Experiment Result.....

Page No.

Known as a primary standard substance and the solution is called a primary standard solution.

SECONDARY STANDARD SOLUTION :

The substance whose solution cannot be prepared directly by weighing its definite amount and then dissolving in definite volume of solvent is called secondary standard substance.

The common secondary standard substances are alkali hydroxide, inorganic acids and KMnO_4 , etc.

NORMALITY (N) :

The normality of a solution is the number of gram-equivalent of the solute per litre of the solution.

$$N = \frac{\text{No. of grams equivalent of solute}}{\text{volume of the solution in } 1000\text{mL}}$$

PROCEDURE :

1. Prepare NaOH solution
2. Rinse and fill the burette with N/5 oxalic acid solution.
3. Pipette out 10mL of above prepared NaOH solution in a conical flask.
4. Add 2-3 drops of phenolphthalein indicator solution

Teacher's Signature : _____

OBSERVATION TABLE:-

S.NO.	VOL. OF NaOH SOLN (vi)	BURETTE READING		VOL OF STD. SOLN OF OXALIC ACID (vi)
		INITIAL	FINAL	
1)	10ml	0.0	9.5	9.5 ml
2)	10ml	10.0	19.4	9.4 ml
3)	10ml	20.0	29.3	9.3ml
4)	10ml	30.0	39.4	9.4ml

Name of Experiment.....

Date: 12/01/2021

Page No.

Experiment No.

Experiment Result

5. add oxalic acid solution from burette drop by drop into it till pink colour disappear.
6. Note down the reading.

OBSERVATION:-wt of the empty weighing tube (w_1) gm.wt of the weighing tube with substance (w_2) gm.wt of the substance ($w_2 - w_1$) gmCALCULATION:-

$$N_1 V_2 \text{ (oxalic acid)} = N_2 V_1 \text{ (NaOH)}$$

$$N_1 V_2 = N_2 V_1$$

$$N_1 \times 10 = 0.1 \times 9.4$$

$$N_1 = \frac{0.1 \times 9.4}{10}$$

$$N_1 = 0.094$$

$$\frac{N_1 \times 9.4}{5} = N_1 \times 10$$

$$N_1 = \frac{0.2 \times 9.4}{10}$$

$$= \frac{1.88}{10}$$

$$N_1 = 0.188$$

RESULT :-

Standard solution of N/10 NaOH solution
is prepared. 0.188

~~0.188~~
~~0.188~~

~~0.188~~
~~0.188~~

Teacher's Signature :

Name of Experiment

Date

Page No.

Experiment No.

Experiment Result.....

PRECAUTION:-

Solution should be made up to desired volume after complete dissolving the solute.


During titration the solution should be stirred thoroughly.

~~Do not take mean of burette readings.~~


Teacher's Signature :

EXPERIMENT-2

REACTION :-

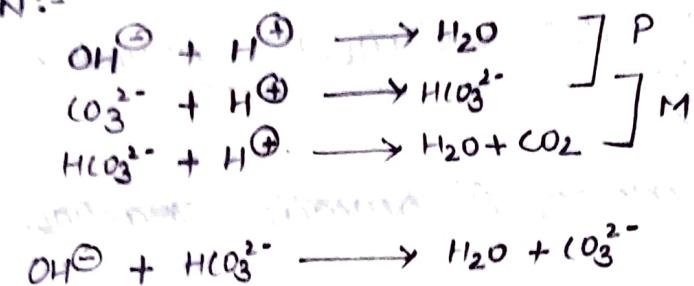

PHENOLPHTHALEIN STRUCTURE :-

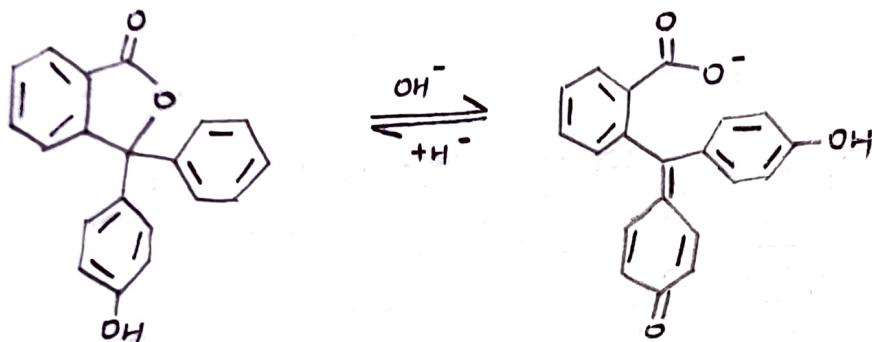

ACIDIC MEDIUM

- 1) COLOURLESS
- 2) BENZENOID STRUCTURE


ALKALINE MEDIUM

- (PINK COLOUR)
- (QUINONOID STRUCTURE)

METHYL ORANGE :-



(BENZENOID STRUCTURE)

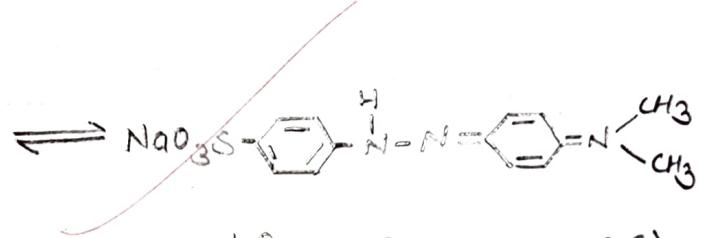

(QUINONOID STRUCTURE)

Date 12/01/2021

Page No.

Name of Experiment

Experiment No.

Experiment Result



EXPERIMENT - 2

OBJECT :- To determine the constituent and amount of alkalinity in the given water sample by titrating it against standard HCl solution ($N/20$) using phenolphthalein and methyl orange as an indicators.

CHEMICALS / APPARATUS :- $N/20$ HCl solution, phenolphthalein and methyl orange indicators, sample solution, burette, pipette, beaker, funnel.

THEORY:- The alkalinity in water is due to the presence of hydroxyl ion (OH^-), carbonate ion (CO_3^{2-}) and bicarbonate ion (HCO_3^-) present in the given sample of water. These can be estimated separately by titrating against standard acid ($N/20$ HCl) using phenolphthalein and methyl orange as indicators.

PROCEDURE:

- 1) Rinse and fill the burette with $N/20$ HCl
- 2) Pipette out 10ml of water sample into a conical flask, add 1-2 drops of phenolphthalein indicator.
- 3) Add HCl from burette drop by drop till the pink colour just disappear.

Teacher's Signature :

$$P = 1 \quad M = 3.1$$

CALCULATION :-

$$\frac{N_1 V_1}{(H_2O)} = \frac{N_2 V_2}{(H_2O)}$$

$$N \times 10 = \frac{1}{20} \times 2P$$

$$N = \frac{1}{200} \times 2 = 0.01$$

FOR CO_3^{2-} -

$$2 \times 1 \times N/20 \times 50 \times 1000$$

$$2 \times 0.05 \times 50 \times 1000$$

$$0.10 \times 5 \times 10^4$$

$$1 \times 10^{-2} \times 5 \times 10^4$$

$$5 \times 10^2 \Rightarrow 500 \text{ mg/L or ppm}$$

For HCO_3^{2-} -

$$2.1 \times 2 \times 0.05 \times 50 \times 1000$$

$$4.2 \times 0.05 \times 5 \times 10^4$$

$$0.21 \times 5 \times 10^4$$

$$1.05 \times 10^4$$

$$10500 \text{ mg/L or ppm.}$$



Name of Experiment.....

Date.....

Experiment No.....

Experiment Result.....

Page No.

- 4) Now add 1-2 drops of methyl orange indicator to above solution. A yellow colour will appear.
- 5) Add $\text{N}/20 \text{ HCl}$ drop by drop to it till the yellow colour just changes to red.
- 6) Repeat this titration until two concordant reading are obtained.

CALCULATION :

 $\text{For } \text{CO}_3^{2-}$

$$N_1 V_1 = N_2 V_2$$

$$(\text{For Water sample}) = (\text{For HCl})$$

$$\text{Volume of HCl used} (V_2) = 2 P$$

$$\text{Strength} = N_1 \times \text{equivalent wt of } \text{CO}_3^{2-} \times 1000 \text{ ppm}$$

$$\text{Alkalinity} = N_1 \times \text{equivalent wt. of } \text{CaCO}_3 \times 1000 \text{ ppm}$$

 $\text{For } \text{HCO}_3^-$

$$N_1 V_1 = N_2 V_2$$

$$(\text{For water sample}) = (\text{For HCl})$$

$$\text{Volume of HCl} (V_2) = (M - P)$$

$$\text{Strength} = N_1 \times \text{equivalent wt of } \text{HCO}_3^{2-} \times 1000 \text{ ppm}$$

$$\text{Alkalinity} = N_1 \times \text{equivalent wt of } \text{CaCO}_3 \times 1000 \text{ ppm}$$

Teacher's Signature :

OBSERVATION TABLE

Sr. No.	VOLUME OF SAMPLE SOLUTN.	BURETTE READING				VOL OF HCl WITH PHENOLPHIT.	VOL OF HCl WITH METHYL ORANGE
		PHENOLPHTHAULEN INITIAL	PHENOLPHTHAULEN FINAL	METHYL ORANGE INITIAL	METHYL ORANGE FINAL		
1)	10	0	1	1	3.1	0.2	1.4
2)	10	0	1.1	1.1	3.0	0.2	1.4

Name of Experiment..... Date.....

Page No. _____

Experiment No..... Experiment Result.....

RESULT:

The alkalinity due to presence of CO_3^{2-} is 500 ppm.

The alkalinity due to presence of HCO_3^{2-} is 10500 ppm

The alkalinity due to presence of OH^- is - ppm.

PRECAUTIONS :

Phenolphthalein indicator should be added first and then methyl orange.

The volume of indicator should same in all the titrations.

The reaction mixture should be shaken properly.

Do not take mean of burette readings.



Name of Experiment.....

Date 19/01/2021

Page No. _____

Experiment No.

Experiment Result.....

EXPERIMENT - 3

OBJECT :- To determine the temporary, permanent and total hardness of water in a given sample of water by titrating it against standard EDTA solution (N/20) using Eriochrome Black-T as an internal indicator.

APPARATUS/CHEMICALS :- Burette, pipette, conical flask, beaker, measuring cylinder, tripod stand, wire gauze, funnel, filter paper, dropper, standard EDTA solution, ammonium buffer solution, Eriochrome Black-T indicator, Hard water.

PRINCIPLE :- When Eriochrome Black-T (indicator) is added to hard water solution at around 10 pH, it gives wine red coloured complex is titrated against EDTA solution. The complex changes from wine red to original blue color showing the end point.

EDTA complex with divalent metal cation ($M = Ca^{2+}, Mg^{2+}$)
The indicator used is a complex organic compound commercially known as Eriochrome Black-T.

Teacher's Signature : _____

SNO.	Vol of hard water	BURETTE READING INITIAL	BURETTE READING FINAL	Vol of EDTA
1)	10ml	0	8.5	8.5 ml
2)	10ml	0	8.5	8.5 ml.

SNO.	Vol of hard water	BURETTE READING INITIAL	BURETTE READING FINAL	Vol of EDTA
1)	10ml	0	4.5 ml	4.5
2)	10ml.	0	4.5ml	4.5

Name of Experiment Date

Experiment No Experiment Result Page No.

PROCEDURE :-

Rinse the burette with N/20 EDTA soln.

Pipette out 10ml of hard water into a conical flask
Add 3ml of basis buffer solution.

Add 2-3 drop of indicator Eriochrome Black-T.
The colour of the soln becomes wine red.

Add drop by drop EDTA soln, wine red to blue colour

Repeat this titration until two concordant reading are obtained.

Take 50ml of hard water in a beaker and boil it until the volume is reduces to nearly half

Filter it and make it to 50ml by adding distilled water.

Pipette out 10ml of this water and titrate until two concordant reading are obtained.

CALCULATION :-

$$N_1 V_1 = N_2 V_2$$
$$(H_2O) \quad (\text{EDTA})$$

$$N \times 10 = \frac{1}{20} \times 8.5$$

Name of Experiment.....
Experiment No.....

Date.....
Experiment Result.....

Page No.

$$N_1 = \frac{0.425}{10}$$

$$n = N_1 = 0.0425$$

$$\therefore \text{Total hardness} = 0.0425$$

$$N_1 V_1 = N_2 V_2$$

$$N_1 \times 10 = 0.05 \times 4.5$$

$$N_1 = \frac{0.2250}{10}$$

$$y = \boxed{N_1 = 0.0225}$$

$$\therefore \text{Permanent hardness} = 0.0225$$

$$\text{Tempranicy hardness} = 0.0425 - 0.0225 \\ = 0.02$$

RESULTS:-

Total hardness = $0.0425 \text{ mg/litre or ppm}$

~~Permanent hardness = 0.0225 mg/litre or ppm~~

Temporary hardness = 0.02 mg/litre or ppm.

Teacher's Signature :

PRECAUTION :-

The glasswares should be properly rinsed with distilled water.

The reaction mixture should be shaken properly.

The end point should be noted correctly.

The pH should be maintained during titration.

The amount of indicator should be same in all titrations.

~~Do not take mean of burette readings.~~

Name of Experiment..... Date..... 02/02/2021

Experiment No..... Experiment Result.....

Page No.

EXPERIMENT-8

OBJECT :- To determine the % of moisture in a given sample of coal by proximate analysis.

MATERIAL REQUIRED :- Coal, Analytical balance, Electric oven, Desiccator, Silica crucible, long legged tong, etc.

THEORY :- Proximate analysis is the simplest type of analysis and it informs about the practical utilization of coal. The analysis since the result have no absolute significance.

MOISTURE coal contain two type of moisture :-

Free or surface moisture :- The moisture which is lost in air drying is called as free moisture or surface moisture.

Internal moisture :- The internal moisture is that type of moisture which is retained by air dried coal. It is rough indicator of the coal maturity

Teacher's Signature : _____

OBSERVATION :-

wt of crucible + dry coal = 15.26 gm

wt of crucible + wet coal = 19.49 gm

wt of crucible + coal sample = 18.98 gm

wt of coal sample = 4.23 gm

% of moisture = 2%

EXPLANATION :- Coal contains water which is removed by heating.

Water is volatile and it volatilizes during heating. It is removed from the sample with the help of heat.

In the above experiment, after heating, the weight of the sample decreased.

Weight loss is due to removal of water.

Water is present in the coal.

OBS

Weig

Weig

weig



OBSERVATION AND CALCULATIONS:-

Weight of empty crucible = n_1 gm = 15.26 gm

Weight of crucible + coal sample = n_2 gm = 19.49 gm

Weight of crucible + coal sample = n_3 gm = 18.98 gm

∴ Weight of coal sample = $(n_2 - n_1)$ gm. = 4.23 gm

∴ Weight of moisture = $(n_2 - n_3)$ gm = 0.51 gm

$$\% \text{ of moisture} = \frac{(n_2 - n_3) \times 100}{(n_2 - n_1)}$$

$$= 0.0425 - 0.0225$$

$$= 0.02.$$

RESULT :-

moisture content in given coal sample

is 2 %.

~~0.02~~
~~0.02~~